
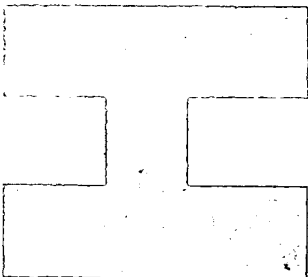
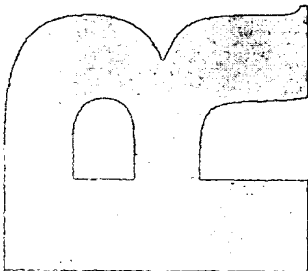
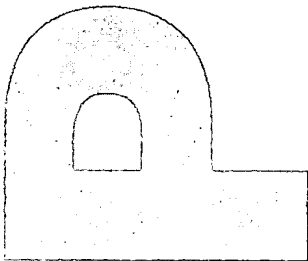
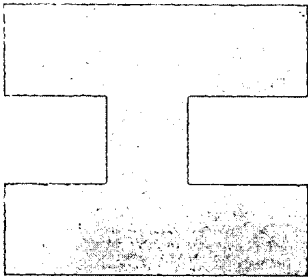


① Prior Art mentioned in specification





P8182US

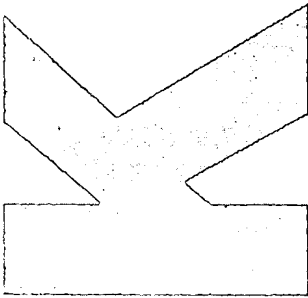

**Anton Paar**

**HR-PHK**  
for NanoSTAR

**Instruction  
Handbook**

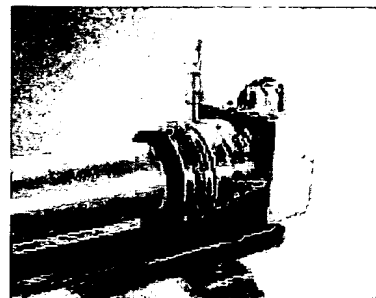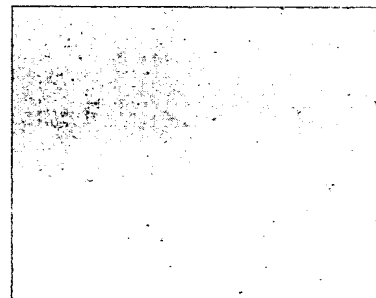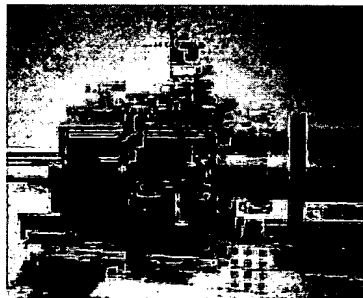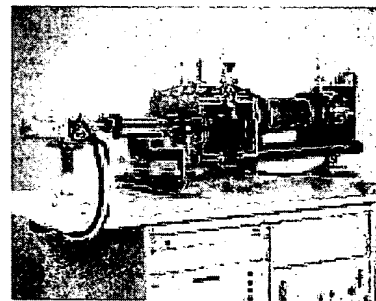

**BEST AVAILABLE COPY**

# **HR-PHK**

## **for NanoSTAR**

Document No.: B23IB06-B

## **Instruction Handbook**

**Anton Paar** GmbH  
Kärntner Straße 322  
A-8054 Graz / AUSTRIA - EUROPE

Phone: +43 316 257-350  
Fax: +43 316 257-257  
E-mail: [doppler@anton-paar.com](mailto:doppler@anton-paar.com)  
Internet: <http://www.anton-paar.com>

Copyright © 1998 **Anton Paar** GmbH, Graz/Austria

While every precaution has been taken in the preparation of this document, **Anton Paar** GmbH assumes no responsibility for errors or omissions.

Neither is any liability assumed for damages resulting from the use of the information contained herein.

**Anton Paar** GmbH does not make a commitment to update the information contained herein.


The instrument may be passed to other users only together with this handbook.

<b>Written by:</b> Thomas A. Binder Richard Görgl	<b>Checked by:</b> Alois Löckner		<b>Released by:</b> Peter Doppler	
	Date:	Signature:	Date:	Signature:

## 5.4 Scattering tube

An ISO-K DN160 flange is used to mount the scattering tube to the sample chamber. The tube is supported on the adjustable stand by two adapter rings. The primary beam stop is fixed at the rear of the scattering tube.


When working in a wide angle range, however, the scattering tube is allowed to be removed. In this case, the primary beam stop has to be fixed directly to the flange connector at the chamber.

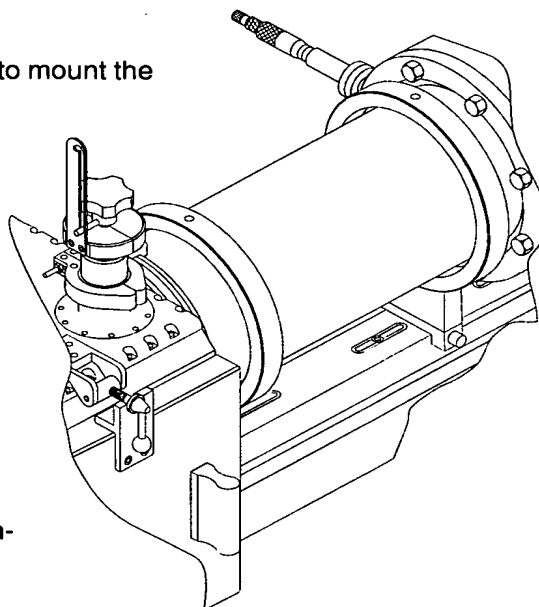

## 5.5 Primary beam stop

The primary beam stop is built into a rigid stainless steel ring. The actual beam stop is made of white gold and its size is matched to the respective diaphragms. A supporting piece is used to attach the stainless steel ring onto the adjustable stand.

The primary beam stop (1) is fixed by two 90° crossed nylon strings and is adjustable over two micrometer drives (2). The primary beam stop is already prealigned – to see the respective values for the micrometer drives refer to the Setting Sticker attached on the sample chamber housing (factory setting).

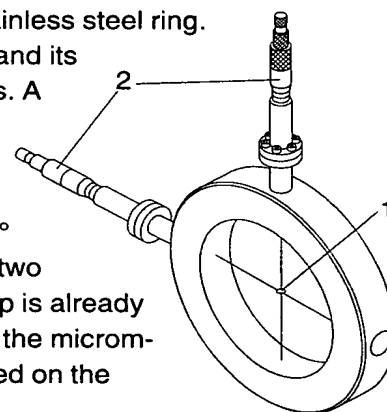

The primary beam stop is flanged to the detector. Two quick-acting couplings (3) are used to press the detector against the O-ring of the primary beam stop ring.

